


CVS[®]
pharmacy

**dental
water jet**

Instruction Manual Model WJ3CVS

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

**BATTERY
OPERATED**
3 AA batteries
(not included)

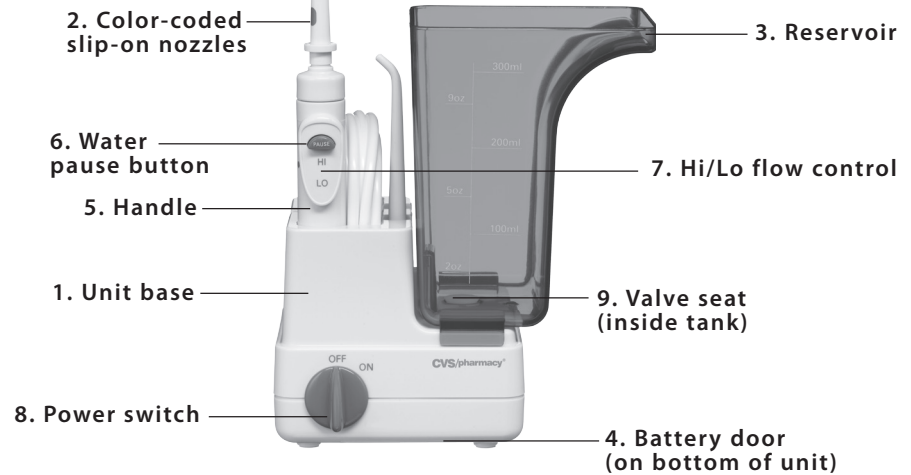

Welcome to a New Level of At-Home Dental Care

The CVS/pharmacy® Dental Water Jet enhances oral hygiene and significantly contributes to healthy gums when used as part of a responsible oral health management program. It is particularly effective at cleaning around orthodontic appliances, fixed and removable bridge work, implants and splints, and it provides valuable gum massage with every use.

Remember: Oral irrigation is a supplement to, not a substitute for, your oral hygiene routine. Ask your dental health care provider to advise you on how the Water Jet can most effectively meet your specific oral hygiene needs.

Initial Preparation

Before using your Dental Water Jet for the first time:

1. Remove battery door (4) and insert 3 AA alkaline batteries. Insert the batteries so that the +/- terminals of the battery match the +/- markings in the battery compartment. Replace door.

NOTE: Use only alkaline batteries. DO NOT MIX old and new batteries. Remove batteries when storing the apparatus for an extended period of time. Clean the battery contacts and also those of the device prior to battery installation. Remove used batteries promptly.

2. Prior to filling the water reservoir, make sure that the water tip is firmly placed onto the handle and into the stand. Once you turn on the Water Jet, it needs several seconds before the water will come out through the tip.
3. Allow full reservoir of lukewarm water to pass through hose, handle and nozzle before first use. (See Instructions for Use, following.)

Instructions for Use

1. Remove handle (5) and slip one color-coded jet tip nozzle (2) onto nozzle support.
2. Remove water reservoir (3) by holding unit base (1) firmly in place with one hand while pulling reservoir (3) straight up with the other.
3. Fill reservoir with water, dental rinse, or mouthwash and reinstall firmly onto the unit.

NOTE: Thoroughly rinse reservoir with clean water after each use with mouthwash or dental rinse, and whenever unit has not been used for an extended period of time.

4. Begin with Hi/Lo flow control at low setting to allow gums to become accustomed to Water Jet cleaning.
5. Press the PAUSE button (6), switch unit to ON, place nozzle tip in mouth, and lean over the sink.
6. Release PAUSE button to begin flow and adjust Hi/Lo control (7) to Hi when comfortable. It usually takes several seconds before the water comes out of the water tip.
7. Hold mouth slightly open and keep elbow at mouth level to let water run back into the sink.
8. Turn handle to direct the jet tip toward tooth surfaces and between teeth to rinse away trapped plaque and particles. The pulsing action of the water massages gums as it cleans the entire mouth.

NOTE: To prevent damage to delicate gum areas, do not point Water Jet directly at gums.

Cleaning and care

If, in the course of use, the valve (9) on the water tank begins to leak, cleaning the valve seat is necessary. To do this, simply use the nozzle tip to push out the valve body, clean the valve seat with warm soapy water, rinse well and replace the valve.

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING

DANGER

- **DO NOT submerge the unit in water or any other liquid.**

WARNING

- Close supervision is necessary when this product is used by, on, or near children or disabled persons.
- Use this product only for its intended use as described in this manual. Do not use attachments unless recommended by the manufacturer.

MEDICAL WARNING

For your safety

- Contact your dentist prior to use if you have had oral or periodontal surgery in the previous twelve months.
- When you start any new oral hygiene program, your gums may initially bleed slightly. Contact your dentist if excessive bleeding occurs or if bleeding continues after two weeks of use.
- Be sure that only the water jet nozzle touches your teeth and gums. Damage to teeth and gums may occur if any other part of the unit touches them.
- This product is designed for oral cleaning only. DO NOT use the unit for any other purpose.
- As with any other new oral hygiene program, individuals with heart problems should check with their physician prior to using the CVS/pharmacy® Dental Water Jet, and should always exercise particular care with any oral hygiene device.
- The CVS/pharmacy® Dental Water Jet is a personal care unit and is intended for household use only.

SAVE THESE INSTRUCTIONS

1 Year Limited Warranty

This product is warranted against defects in materials and workmanship for a period of one (1) year from the date of original purchase or receipt as a gift. This warranty does not cover institutional or commercial usage, or damages caused by misuse, neglect or accident.

Should this product become defective during the warranty period, return it to the place of purchase for repair or replacement, at the retailer's option, without charge to you.

NEITHER THE RETAILER NOR THE MANUFACTURER, IN ANY EVENT, SHALL BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitations may not apply to you. You may also have other rights which vary from state to state.

FOR REPLACEMENT, PLEASE BRING THIS PRODUCT DIRECTLY TO THE PLACE OF PURCHASE WITH THE SALES RECEIPT (PLEASE DO NOT MAIL).